OPERATING INSTRUCTIONS FOR 2x4 50ml & 25mL Tube SpinVessel® VP 418SV2-2x4-50-CC

US Patent #11,623,188
European (BE, FR, DE, NL, CH, SE, & UK) Patent #3887049



SpinVessel® System, with Computer Control feature, for 25mL and 50mL SpinVessel® Tubes:

SpinVessel® System VP 418SV2-2x4-50-CC includes a Motor Unit with Spin Bases, as well as a Controller and Power Supply (not shown). Shown with VP 830SV-50 SpinVessel® Tubes, which are sold separately.

Table of Contents

Components ---------- pg. 3


Setup ---------- pg. 4


Operation ---------- pg. 5

Product Maintenance ---------- pg. 8

Warranty ---------- pg. 9

Components

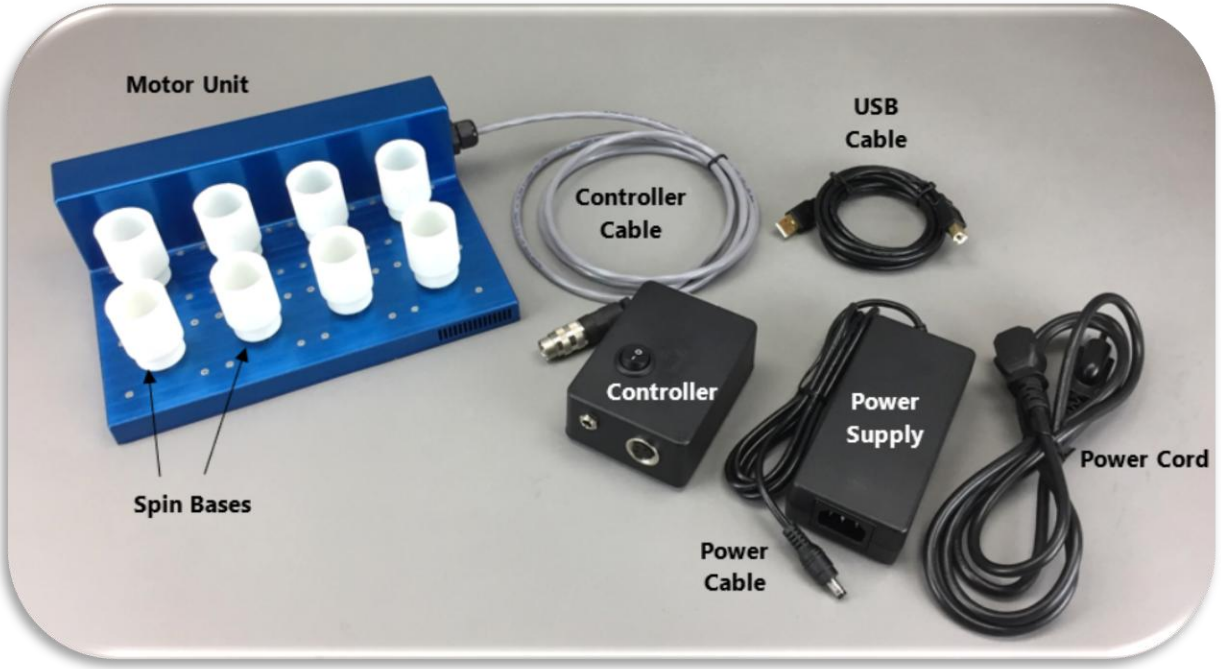
Motor Unit
USB
Cable
Controller
Cable
Controller
Power
Supply
Power Cord
Power
Cable
Spin Bases

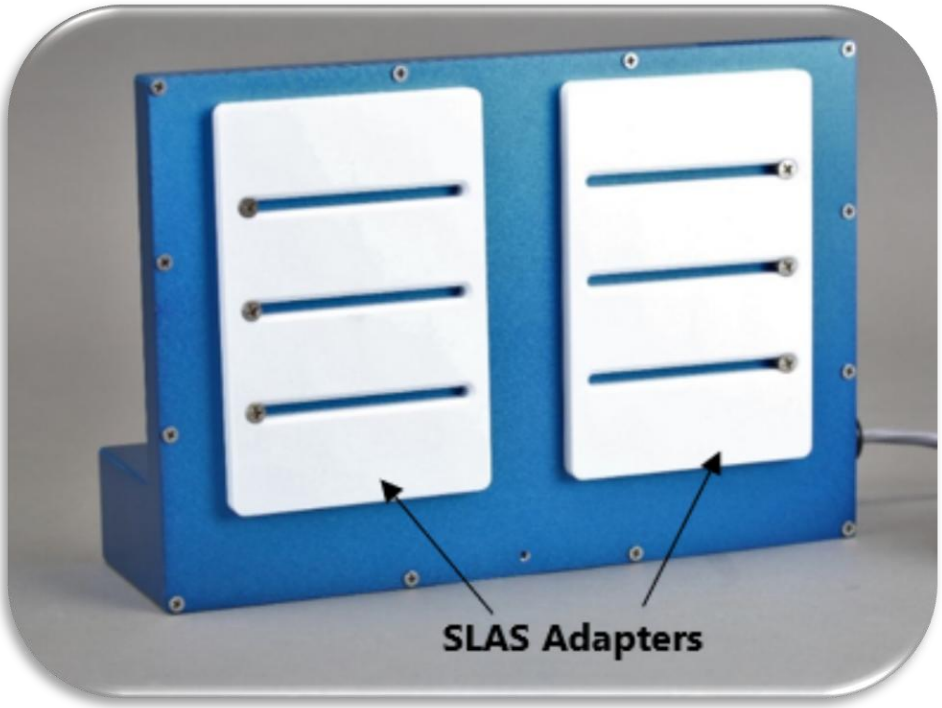
SLAS Adapters

POSITION THE SpinVessel® SYSTEM ON A LIQUID HANDLER

1. Place Motor Unit on a sturdy lab bench, table, or robotic workstation.

    a) SLAS dimensioned adapters are part of the system and enable users to position the SpinVessel® system on a robotic liquid handler.
    b) Adapters have adjustment slots to allow for different spacing of microplate positions on various liquid handling workstations.


Warning: Do not use the Spin Cups as a handle to move the SpinVessel. Hold the Motor Unit with both hands to position the SpinVessel where desired.

2. Since the controller cable is 6 feet long, the controller and computer can be placed at a distance from the Motor Unit.

3. Connect the Power Supply to the Controller Box

4. Connect the Controller Cable from the motor base to the Controller Box

5. Connect the USB Cable to the Controller Box and also to the PC you'd like to use to issue commands

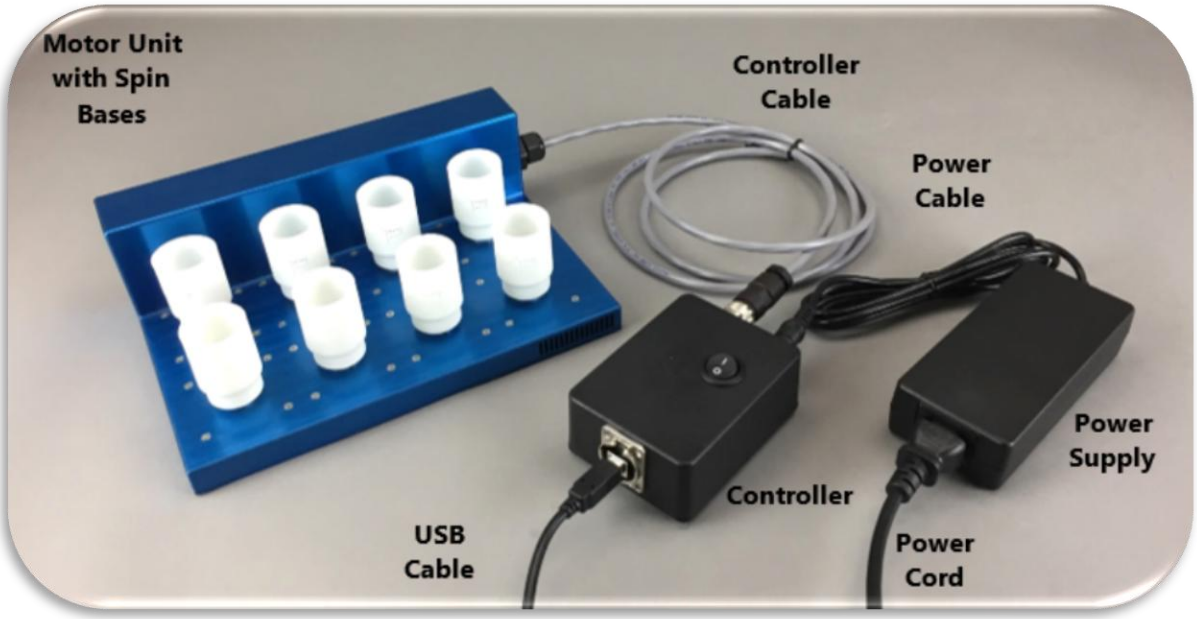

POWERING UP THE SpinVessel® SYSTEM

1. The Power Switch (I/O) is located on the top of the Controller.
2. To power up the Controller, depress the "I" of the switch.

OPERATING THE SpinVessel® SYSTEM WITH COMPUTER CONTROL


The SpinVessel® Models with "-CC" in the part number are controlled using a computer. In this configuration, the control unit does not have a manual Touch Screen Controller, however, V&P Scientific does offer a manual touch screen controller option (VP 418SV-MC) sold separately. The Controller for the computer-controllable SpinVessel® is connected by a USB communication cable to a computer. The Controller accepts simple ASCII commands to set the speed of the rotations, the number of rotations (including fractions of a rotation), the duration of pause before reversing direction and the degree of ramping to set speed.

Software details

1. The SpinVessel® Controller uses an internal CH340G USB to Serial Converter (set to 9600,8, N,1) which will require driver software. Windows may already have the needed USB driver software installed. If not, download the "CH340 USB DRIVER" from the provided USB flash drive and run "CH341SER.exe" to install it (Windows typically auto-installs once downloaded). If your company policy does not allow the use of the provided USB stick, this driver is also available for download from the internet (contact your IT department for assistance if necessary).
2. The basic command language for this computer-to-stirrer communication is ASCII.
3. Commands can be sent from the computer using a standard terminal program like HyperTerminal or Realterm. Or, if needed, download the V&P Scientific Serial Terminal program from the provided USB flash drive. If your company policy does not allow the use of the provided USB stick, we suggest asking your IT department for a recommendation.
4. Alternatively, the ASCII commands can be sent by an automated liquid handler's software to control the SpinVessel®. Since there are many different automated liquid handlers available, please contact a technical representative from the company that makes the liquid handler.

Connecting and sending commands

1. Set up the system as described in the Set-Up section of the SpinVessel® Technote with the USB cable connected to the computer to be used.
2. Open the Serial Terminal program (V&P Scientific version or other).
3. Turn on the SpinVessel® instrument as described.
4. Select the appropriate COM port and click CONNECT to establish a line of communication between the instrument and the computer.
5. Create ASCII command strings as outlined in the diagram shown below to control the instrument with the following parameters:
    a) Speed of rotation in revolutions per minute (RPM). Do not exceed 650 RPM.
    b) Number of rotations before reversing direction.
    c) Fraction of a rotation (input 0-9 for fractions of a rotation).
    d) Pause, in milliseconds, before reversing direction.
6. Hit ENTER to send command and start rotation of the SpinVessel® instrument.
7. To stop, type <0> and hit ENTER.

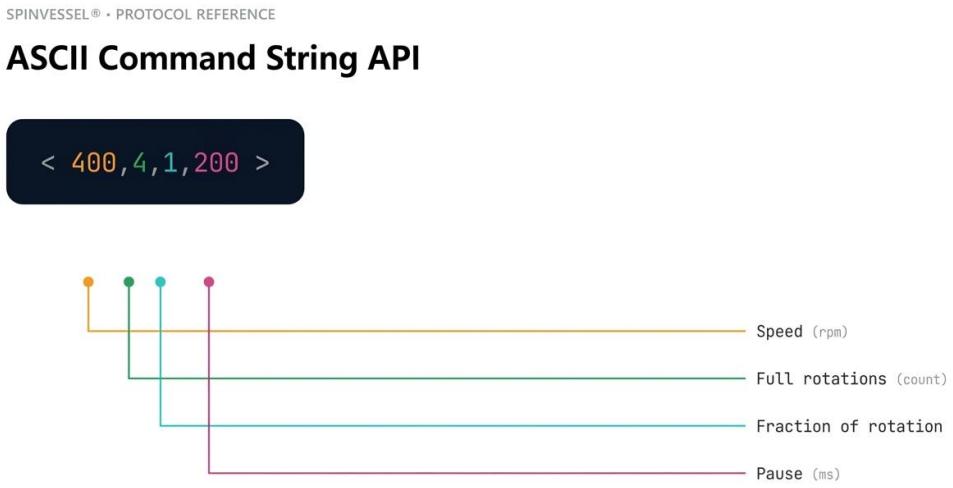

USING the 25mL or 50mL SpinVessel® SYSTEM


The optimal stirring mode is dependent upon the application and needs to be empirically determined. Factors to consider in determining optimal 25mL or 50mL SpinVessel operation are the density of the particulates, their size and shape, as well as the volume and viscosity of the liquid. Start with 400 RPM and 4 rotations before reversing direction, then test to determine the appropriate RPMs and number of rotations for the particulate's density, fragility, and solution viscosity. We have found that smaller circumference SpinVessels® require higher RPMs than larger circumference SpinVessels, because the critical factor is the linear speed generated at the outer edge of the vessel. See page 501 of our SLAS Technology paper for a full discussion of this phenomenon.
The following settings are ones to use when starting out with the VP 418SV2-2x4-50-CC:

Speed: 400 RPM Rotate: 4 times Pause: 200 milliseconds

It is always a good idea to practice with water in the SpinVessel® to determine safe RPM, rotation, and pause settings before using expensive reagents.

The VP 418SV2-2x4-50-CC SpinVessel® is designed to be used with VP 830SV-50 or VP 830SV-25 SpinVessels® tubes. Please contact V&P Scientific for pricing.

See below for an example of low versus high RPM for a solution of yeast cells in a large-circumference, 6-liter SpinVessel®. Please also note that the required speed (RPM) for optimal mixing is directly related to the diameter of the SpinVessel tube. The smaller the diameter, the faster the RPMs need to be to effect adequate mixing. Note that, when given a longer mixing time, even the slower speed (20 RPM for gentler mixing) was able to suspend the yeast cells.

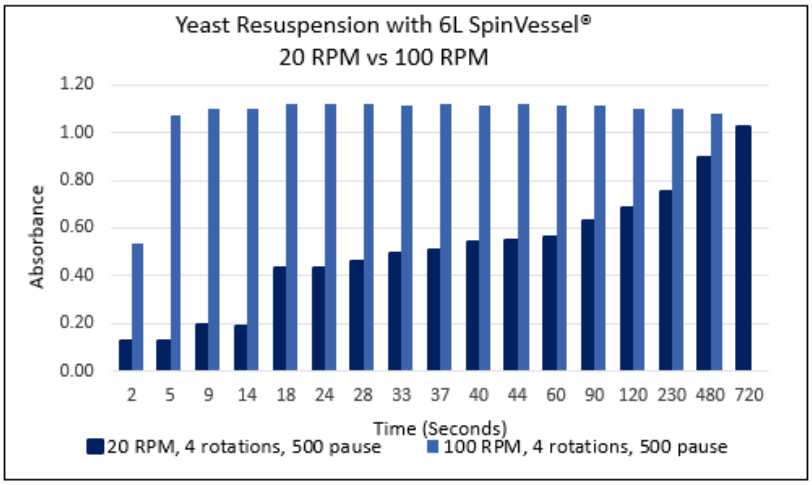

GENERAL PRODUCT CARE

When not in use, turn the power switch off.
Do not place the Controller in chambers with temperatures above 40°C.
To clean the SpinVessel, wipe down with a cloth and mild detergent, followed by a water wipe.
Do not immerse SpinVessel® Motor Unit or Controller in liquid.

SAFETY PRECAUTIONS

The use of motor controls, like that of all utilization of concentrated power, is potentially hazardous. The degree of hazard can be greatly reduced by proper design, selection, installation, and use, but all hazards cannot be completely eliminated.

The following safety precautions must be observed during all phases of installation, operation, service, and repair of this motor control product. Failure to comply with these precautions or with specific warnings elsewhere in this manual violates safety standards of design, manufacture and intended use of the products. V&P Scientific assumes no liability for the customer's failure to comply with safety requirements and practices.

Never reach into the consumable vessel by hand while the unit is running. Only precisely controlled robotic applications have the necessary accuracy to avoid spinning projections. Only do so after the movement has come to a complete stop.

Secure all loose clothing, hair, and jewelry before operating the equipment. Loose items can easily become caught in the rotating mechanism, leading to entanglement injuries.

Disconnect the power supply before performing any maintenance, cleaning, or adjustments. Unplug electric mixers or lock out the power source to prevent the machine from accidentally starting while you are working on it.

WARRANTY

V&P Scientific, Inc. warrants this product to be free from defects in material and workmanship when used under normal laboratory conditions for one year. This warranty begins from the date of delivery by V&P Scientific.

In the event this product fails under normal laboratory conditions within the specified period of time because of a defect in material or workmanship, V&P Scientific will, at its option, repair or replace the product. Damage to the product caused by user negligence is not covered.

This warranty is made in lieu of other warranties expressed or implied including the warranties of merchantability and fitness for a particular purpose. V&P Scientific shall not be liable for loss or damages arising from the use of these products nor for consequential damages of any kind.

Please keep the special shipping carton in case the unit needs to be shipped back to V&P Scientific. Contact V&P Scientific for return authorization and shipping instructions or for any other assistance at 858-455-0643 or sales@vp-sci.com.